## Mod. 80

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ENGLISH

## ATTENTION

- Read carefully these instructions before using the appliance for the first time and keep it for future enquiries.
- Keep all plastic bags and packaging components out of the reach of children. They are potentially dangerous.
- Do not immerse the appliance in water nor place it in wet surfaces.
- Do not maintain exposed to extreme temperatures (5 - 35 °C).
- Firmly grasp the item.
- Clean with a slightly damp cloth. Do not use abrasive or chemical products.

## INSTRUCTIONS OF USE

## BATTERY INSERTION

- The scale works with 3 batteries of 1,5 V. type **AAA** (Not included).
- Open the battery housing lid.
- Insert the battery and make sure the poles are correctly positioned and then close the battery door again.
- Remove the batteries if you are not going to use the appliance for a long time.

## MEASURING

- Press **"SET"**, the appliance will turn on automatically.
- You can choose the measuring system in centimetres **"cm"** or feet **"ft"** by keeping pressed **"SET"**. Press **"SET"** once to change the measuring. Once chosen and after 3 seconds the display will show the main screen.

• **Method 1: Child height measuring:**

- 1-Place the child to stand upright against the wall and place the appliance over the child's head as showed on **Figure 1**. Once correctly placed the child can move it away.
- 2-Keeping the appliance hold in the wall, press **"SET"** button till listening a beep. The measurement will be showed in the display.

• **Method 2: Measuring yourself:**

- 1-Keep the lower part of the appliance hold in the wall with your hand.
- 2-Place yourself straight against the wall and place the appliance on your head in horizontal position.
- 3-Do not move it when going back and press the button, (**Figure 2**) the result will be showed in the display after listening a buzzing.

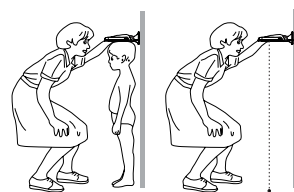
FIGURE 1

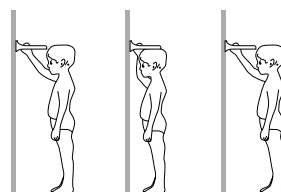
FIGURE 2

## YOU MUST BEAR IN MIND

- **IMPORTANT:** to ensure the accuracy, there can not be a prominent barrier including baseboards within 80 cm diameter range. Please check if there is any edge and move into a wide area when the result is not correct.
- Avoid exceeding the maximum measurement of this item (**200 cm**). If the maximum measurement is exceeded the message “Err” will appear on the display.
- You must proceed to change the battery when the “Lo” symbol is showed on the display. To do that, open the cover placed in the back side of the appliance. Change the batteries.
- The item is a sensible electronic instrument. Some appliances of high frequency as mobiles, radio stations, remote controls and microwaves, can cause interferences if the scale is used nearby of them. Keep it far away of these points if you notice irregularities on the display.
- It is not apt for commercial transactions or medical purposes.
- Designed for domestic use only.

## ENVIRONMENT PROTECTION

- Dispose the appliance according to the local garbage regulations.
- Dispose the batteries according to the local garbage regulations.
- Never throw away the appliance or the batteries to the rubbish. So you will help to the improvement of the environment.
- The batteries must not be mixed with other domestic residues.

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**DECA,S.A.**